# beurer

## GS 80

**(D) Glaswaage**  
Gebrauchsanleitung

**(GB) Glass scale**  
Instruction for use

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## Contents

|                                          |   |                                              |    |
|------------------------------------------|---|----------------------------------------------|----|
| 1. Getting to know your instrument ..... | 7 | 7. Measuring your weight .....               | 9  |
| 2. Safety notes .....                    | 7 | 8. Cleaning and maintaining the device ..... | 10 |
| 3. Unit description .....                | 8 | 9. Disposal .....                            | 10 |
| 4. Initial use .....                     | 8 | 10. What if there are problems? .....        | 10 |
| 5. Settings .....                        | 8 | 11. Technical specifications .....           | 11 |
| 6. Changing the batteries .....          | 9 |                                              |    |

## Included in delivery

- Glass scale
- Outdoor sensor
- 4 x AA batteries (scale)
- 2 x AAA batteries (outdoor sensor)
- These instructions for use

## Dear Customer,

Thank you for choosing one of our products. Our name stands for high-quality, thoroughly tested products for the applications in the areas of heat, weight, blood pressure, body temperature, pulse, gentle therapy, massage and air.

Please read these instructions for use carefully and keep them for later use, be sure to make them accessible to other users and observe the information they contain.

Best regards,

Your Beurer Team

## 1. Getting to know your instrument

### Functions of the device

This scale allows you to weigh yourself and also displays the time, inside and outside temperature as well as a weather forecast.

The scale includes the following features:

- Body weight measurement
- Time
- Inside and outside temperature
- Weather forecast
- Alarm function
- Switch between kilograms (kg), pounds (lb) and stones (st)
- Automatic switch-off function
- Low battery indicator

## 2. Safety notes

### Warning:

- Never step onto the outer edge of the scale with one foot as this could cause the scale to tip!
- Keep packaging material away from children (risk of suffocation).
- Swallowing batteries can be extremely dangerous. Keep the batteries and scale out of the reach of small children. Should a battery be swallowed, seek medical assistance immediately.
- Batteries must not be recharged, taken apart, thrown into an open fire or short circuited.
- Batteries can contain toxins that are harmful to health and the environment. Always dispose of batteries in accordance with applicable legal regulations. Do not dispose of batteries with the normal household waste.


## General notes

- The device is suitable only for self-testing, not for medical or commercial purposes.
- The scale's maximum capacity is 150 kg (330 lb/24 st).
- Check the batteries before submitting a repair query and replace them if necessary.
- Repairs must only be carried out by Customer Services or authorised suppliers.
- If you should have any questions about using our devices, please contact either your retailer or Customer Services.


## 3. Unit description

- 1 Weight display
- 2 Time
- 3 Weather forecast
- 4 Inside temperature
- 5 Outside temperature
- 6 Channel
- 7 Outdoor sensor
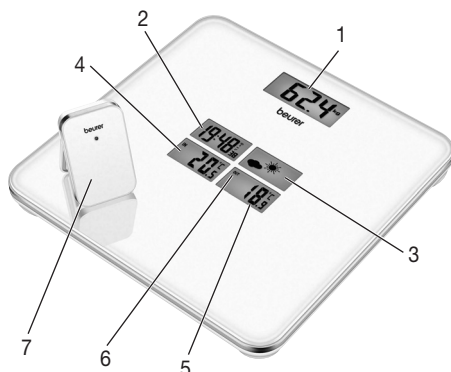
## Battery compartment

- 1 MODE button
- 2 CHANNEL button
- 3 ALARM button
- 4 SNZ button
- 5 ▲ button
- 6 ▼ button
- 7 Reset button
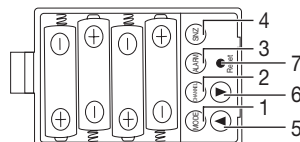
## 4. Initial use

Remove the batteries from the protective packaging and insert them into the scale and then into the outdoor sensor within 60 seconds, ensuring that the battery polarity is correct in all cases.

After inserting the batteries, "--.--" appears in the outside temperature LCD field and the radio symbol flashes. The scale now receives the signal from the outside temperature sensor. The transmission is complete when the outside temperature is displayed.

**If the scale does not function, remove the batteries completely and reinsert them. Ensure that the correct battery polarity is observed. You can also press the "Reset" button. Press the CHANNEL button to select the first channel (CH1), as the outdoor sensor always searches for the first channel.**

## 5. Settings

### Setting the time

In standby mode (time displayed), press the MODE button for 3 seconds. Now use ▲ and ▼ to choose between 24- and 12-hour mode. Press the MODE button to confirm your selection. Now set the time (hours and minutes) using ▲ and ▼, in each case confirming your entry with the MODE button.

### Setting the alarm clock

Press the MODE button in standby mode. "AL" appears in the display (2). Press and hold the MODE button for 3 seconds. You can now set the desired alarm time using ▲ and ▼. Press the MODE button to confirm your entry.

**Activating/deactivating the alarm clock**



Press the ALARM button in standby mode. The alarm symbol  appears in the display (2). Pressing the ALARM button again deactivates the alarm clock. An accelerating acoustic signal sounds at the set alarm time. Press the SNZ button to activate the snooze function; in this case, the alarm sounds again after 5 minutes. Press any button to deactivate the alarm clock.

**Setting the temperature unit**

In standby mode, press the  button to select between °C and °F.

**Configuring the weather station**

Start by setting the current weather situation. Press and hold the  button until an acoustic signal sounds. Now use the  button to select the weather situation 12 hours previously. Confirm using the  button. The station displays the next change of weather in the next 12 hours. The scale uses a range of different weather symbols to provide weather information in as much detail as possible. Depending on the weather forecast, the scale displays a sequence of flashing weather symbols:
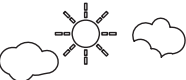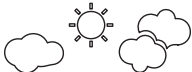
|                 | 1.                                                                                 | 2.                                                                                 | 3.                                                                                 | 4.                                                                                 |
|-----------------|------------------------------------------------------------------------------------|------------------------------------------------------------------------------------|------------------------------------------------------------------------------------|------------------------------------------------------------------------------------|
| Sunny           |   |   |   |   |
| Partly Cloudy   |   |   |   |  |
| Cloudy          |   |   |   |  |
| Rainy           |   |   |   |  |
| Torrential Rain |   |   |   |  |
| Snowy           |  |  |  |  |

 The weather forecast is based on changes in atmospheric pressure. Pressure values are compared on an hourly basis. If the atmospheric pressure rises, an area of high pressure can be expected. The indicated forecast changes towards the “sunny” end of the weather scale. The opposite applies in the case of declining atmospheric pressure.

**6. Changing the batteries**

The scale is equipped with a low battery indicator. If the batteries become too weak, “Lo” appears in the display (1) and the scale switches off automatically. In this case, the batteries must be replaced as soon as possible (4 x 1.5 V, type AA). The batteries in the outdoor sensor must be replaced when the scale no longer receives a radio signal (2 x 1.5 V, type AAA). Press and hold the CHANNEL button for 3 seconds to search for the outdoor sensor. Ensure that channel 1 (CH1) is set on the scale.

**7. Measuring your weight**

Place the scale on an even, hard surface (no carpet); a hard surface is crucial to accurate measurements. Tap the scale with your foot; “RL” briefly appears in the display, followed by “0.0 kg”. Step onto the scale. Stand still with your weight distributed equally on both feet. The measurement is displayed shortly afterwards.

If you step off the standing surface, the scale switches off after 10 seconds. The measurement remains visible during this period.

If the maximum capacity of the scale (150 kg) is exceeded, “0-Ld” appears in the display.

## 8. Cleaning and maintaining the device

The device should be cleaned from time to time.

This should be done using a damp cloth and, if necessary, a small amount of detergent. Do not use any abrasive cleaning products and never submerge the device in water.

### Important:

- Never use abrasive solvents or cleaning agents!
- Never submerge the device in water!
- Do not clean the device in the dishwasher!

## 9. Disposal

Batteries must not be disposed of with household waste. As a consumer, you are required by law to recycle used batteries. You can recycle your old batteries at public collection points in your community or wherever batteries of the relevant type are sold.

### Note:

The codes below are printed on batteries containing harmful substances:

Pb = Battery contains lead,

Cd = Battery contains cadmium,

Hg = Battery contains mercury.

Observe the local regulations for material disposal.

Dispose of the unit in accordance with EC Directive **2002/96/EC – WEEE** (Waste Electrical and Electronic Equipment). If you have any questions, please contact the local authorities responsible for waste disposal.

## 10. What if there are problems?

If the scale encounters an error during measurement, the following is displayed:

| Display | Cause                                                       | Solution                                                 |
|---------|-------------------------------------------------------------|----------------------------------------------------------|
| 0_Ld    | Maximum capacity of 150 kg exceeded.                        | Loads must be limited to 150 kg.                         |
| Lo      | Batteries are almost empty.                                 | Replace the batteries in the scale.                      |
| CAL     | The scale is finding its zero point.                        | Wait until “0.0 kg” is displayed.                        |
| LL:L    | Inside temperature < -10 °C<br>Outside temperature < -50 °C | Inside or outside temperature outside measurement range. |
| HH:H    | Inside temperature > 50 °C<br>Outside temperature > 70 °C   | Inside or outside temperature outside measurement range. |

### No measurement possible

| Possible faults                                                                                                                                          | Solution                                                                                                                                                       |
|----------------------------------------------------------------------------------------------------------------------------------------------------------|----------------------------------------------------------------------------------------------------------------------------------------------------------------|
| The scale was not activated before stepping onto it. If you stand on the scale before the display shows “0.0 kg”, the scale will not function correctly. | Correctly activate the scale (wait for “0.0 kg” to be displayed) and repeat the measurement.<br>If necessary, wait 30 seconds and then repeat the measurement. |

| Possible faults                                                | Solution                                                                                                                                                                                                                                                                                                        |
|----------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Radio connection between scale and outdoor sensor interrupted. | Reinsert the batteries in the scale and outdoor sensor; the scale emits an acoustic signal when the batteries are inserted.<br>Ensure a sufficient distance from the following electromagnetic sources: electricity pylons, TV or computer radio systems, fitness devices, mobile telephones, weather stations. |
| The outdoor sensor is too far from the scale.                  | Reduce the distance between the scale and outdoor sensor to a maximum of 30 meters (open field).                                                                                                                                                                                                                |
| The batteries in the scale are empty.                          | Replace the batteries in the scale (see page 9).                                                                                                                                                                                                                                                                |
| The batteries in the outdoor sensor are empty.                 | No weather forecast is displayed. Replace the batteries in the outdoor sensor.                                                                                                                                                                                                                                  |

We hereby guarantee that this product complies with the European R&TTE Directive 1999/5/EC.  
Please contact the specified service address to obtain further information, such as the CE Declaration of Conformity.

## 11. Technical specifications

|                        |                                |
|------------------------|--------------------------------|
| Measurement range      | 2–150 kg, 100 g graduation     |
| Inside temperature     | -10 °C – +50 °C                |
| Outside temperature    | -50 °C – +70 °C                |
| Temperature graduation | 0.1 °C                         |
| Radio transmission     | 433 MHz                        |
| Range                  | 30 metres (open field)         |
| Inside temperature     | 30 second measurement interval |
| Outside temperature    | 57 second measurement interval |

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